NDODONTICS

Create it



iPexII

Advanced Apex Locator

Using the results of numerous clinical trials and verification tests as a base, NSK created SmartLogic, a next-generation Apex Locator compatible with diverse tooth forms, and that provides superior accuracy in detecting the root apex.

At the biologically crucial root apex - the Zone of Precision - the iPexII is able to provide a clear, accurate image of the file tip's location on its LCD panel.

Stylish & Compact Body

Lightweight and with a base just 6 cm in diameter, the iPexII can quickly be set up and used in virtually any environment. The LCD panel is pitched at the ideal angle that allows clinicians to easily and without distraction monitor conditions ensuring safe and smooth root canal procedures.

0.0 NSK iPexII

Next-Generation SmartLogic Controller

NSK's SmartLogic controller combines two advanced technologies to accurately analyse signals. iPexII SmartLogic technology removes virtually all signal interference from the root canal itself to precisely detect the location of the apex.


- iPexII SmartLogic alternates between two separate frequencies to take measurements. By reducing the noise that interferes with signal analysis, it generates a high-precision wave system that does not require a filtering circuit.
- iPexII SmartLogic uses the whole wave system to analyse the extracted frequency. This produces steady measurement analysis even under varying conditions in the root canal while retaining a constant link to file insertion progress.

The 3 colour LCD visual monitor screen

▶ Measurement Begins





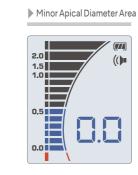

Overshooting the Apex



Monitoring the Canal Depth

For integration into Endodontic instruments (ENDO-MATE TC2 / ENDO-MATE AT / Contra-angle) The NSK endodontics products and iPexII ensure safe, precise and efficient root canal procedures.



iPexII Complete Set

iPexII Complete Set	Y1002208
MODEL	ORDER CODE

Contents
• Control unit • iPexII Probe (1.8 m) • iPexII File Clip (3 pcs.)

• Lip Hook (3 pcs.) • AAA manganese dry batteries (3 pcs.) • iPexII Tester

- Specifications
 Power Supply: for 1.5 V AAA (LR03 or R03) batteries
- Dimensions : W 60 x D 60 x H 86 mm
- Weight: 76 g (unit without batteries)



Spare Parts

ENDO-MATE TC2

MODEL	PRODUCTS	ORDERCODE	MODEL	PRODUCTS	ORDERCODE	MODEL	PRODUCTS	ORDERCODE
ENDO-MATE TC2 Motor	S (5/2/0)	E346	Battery Charger		U614	Cord for MPA / MPAS		C876021
Ni-MH Rechargea Battery (Pack of 2)	ble 🚺	Y900103	Power Cord	6 6	U439551	Connection Cord Clip		C876509

^{*}Use only a rechargeable battery, size AAA Ni-MH (Nickel Metal Hydride). Never attempt to use a non rechargeable battery.

ENDO-MATE AT

MODEL	PRODUCTS	ORDER CODE	MODEL	PRODUCTS	ORDER CODE
ENDO-MATE DT Adaptor		U480	Cord for MPA / MPAS		C876021
Handpiece Stand		Z095205	Connection Cord Clip		C876509

ENDO-MATE TC2 / ENDO-MATE AT Heads

MODEL	PRODUCTS	ORDERCODE	MODEL	PRODUCTS	ORDERCODE
MP-Y Miniature Head, for Ni-Ti files (ø2.35), Push Button Chuck		C547	MPA-Y Connection for Apex locator (for long shank), Miniature Head, for Ni-Ti files (82.35), Push Button Chuck		C548
iPexII					
MODEL	PRODUCTS	ORDER CODE	MODEL	PRODUCTS	ORDER CODE

13

MODEL	PRODUCTS	ORDER CODE	MODEL	PRODUCTS	ORDER CODE
PexII File Clip		U1109351	iPexII Probe		U1109352
PexII Tester	1	U1109353	Lip Hook		U501513